

TAKTVOLL

Intelligent Manufacturing in China



Vessel Sealing Generator

For Veterinary Use


ES-100VL

Both affordable and portable, easy-to-use. The platform is intended for open and laparoscopic surgical procedures.

Vessel Sealing Generator

ES-100VL For Veterinary Use

- Multi-mode with Ligasure.
- Automatic shutdown when sealing is done, without operator' s involvement.
- Intelligent power regulation, with lower energy to achieve best sealing on target tissue.

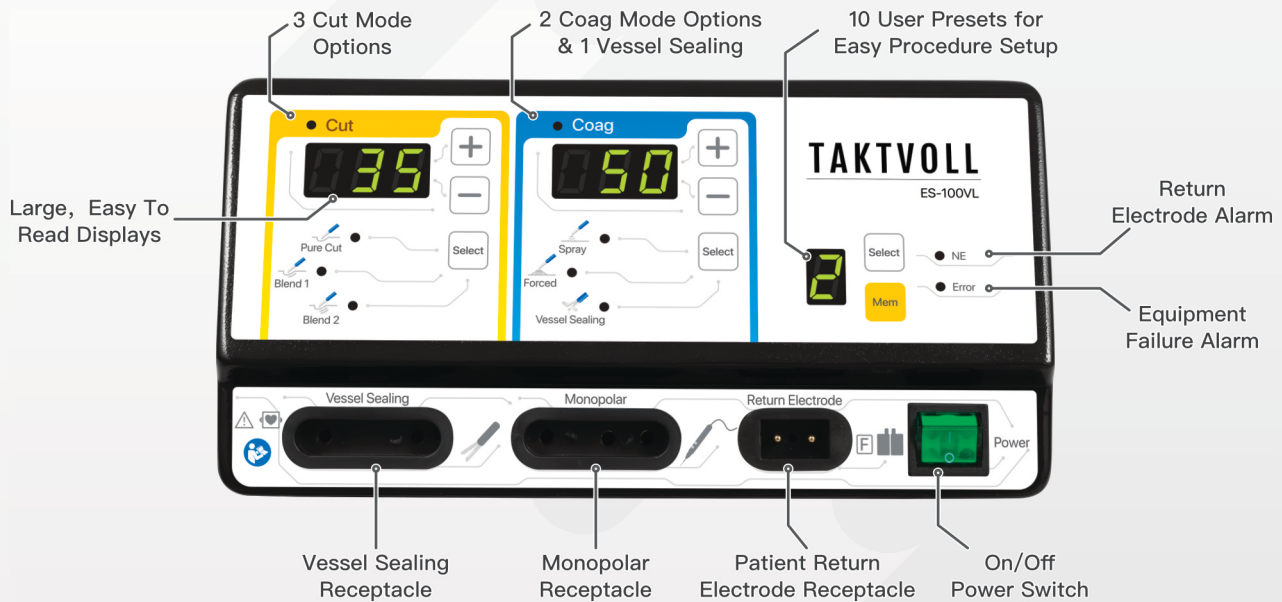

ES-100VL Feature

3 Cut Mode

PURE

Provides a standard Cutting effect.

BLEND 1

Provides a combination of 70% Cutting and 30% Coagulation

BLEND 2

Provides a combination of 50% Cutting and 50% Coagulation

2 Coag Mode

SPRAY

Designed to provide hemostasis in a broader area with less depth of tissue damage when compared to desiccation.

FORCED

Fast and effective standard coagulation with moderate thermal penetration. Slight carbonization effects.

1 Vessel Sealing

Vessel Sealing

Seal and cut vessels up to and including 7mm in diameter.

Sealing vessels and tissue bundles without dissection or isolation , easy to operate.

TAKTVOLL®

Beijing Taktvoll Technology Co., Ltd.

Address: Room 10a-302, No. 13, South Fourth Street, Jingsheng, Jinqiao Science And Technology Industry Base, Tongzhou District, Beijing, China.

Tel: +86(10) 56495579 / 56495580 Web: www.taktvoll.com.cn E-mail: Support@taktvoll.com.cn