

TAKTVOLL®

**NEW PRODUCT!** ES-100V PRO LCD  
Touchscreen Electrosurgical System  
**ES-100V Pro**




The ES-100V Pro satisfies the veterinarian's demands with precision, safety, and reliability.



# Vessel Sealing Generator

## ES-100VL Pro For Veterinary Use

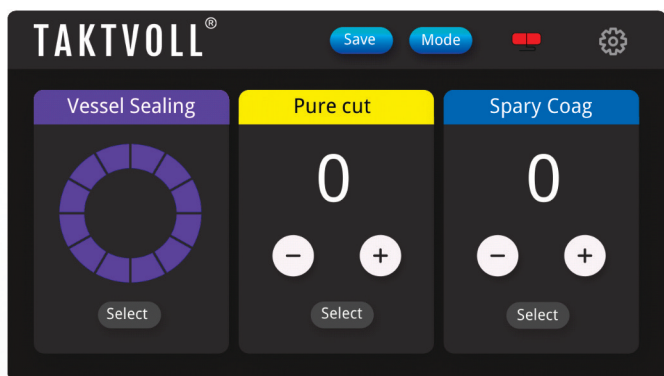

- Glass full version design, which is extremely tense and beautiful
- Color touch screen operation panel, flexible and easy to operate
- 7 Working modes
- Large blood vessel sealing function
- Electrosurgical Pens and Foot Switch Control
- Recently Mode, Power and Other Parameters
- Quick Adjustment of Power
- Cut and Coag in an Intermittent Manner



### ES-100V Pro Feature

| 3 Cut Mode     |   | 2 Coag Mode   |   | 1 Vessel Sealing  |  |
|----------------|---|---------------|---|---|--|
| <b>PURE</b>    | Provides a standard Cutting effect.                       | <b>SPRAY</b>  | Designed to provide hemostasis in a broader area with less depth of tissue damage when compared to desiccation. | <b>Vessel Sealing</b>   |  |
| <b>BLEND 1</b> | Provides a combination of 70% Cutting and 30% Coagulation | <b>FORCED</b> | Fast and effective standard coagulation with moderate thermal penetration. Slight carbonization effects.        | Seal and cut vessels up to and including 7mm in diameter. Sealing vessels and tissue bundles without dissection or isolation , easy to operate. |  |
| <b>BLEND 2</b> | Provides a combination of 50% Cutting and 50% Coagulation |               |   | <b>1 Bipolar</b>  |  |
|                |   |               |   | <b>Standard</b>   |  |
|                |   |               |   | Bipolar coagulation which gently coagulates tissue for minimal tissue damage reducing carbonization and sticking.                               |  |