


# TAKTVOLL®

Intelligent Manufacturing in China

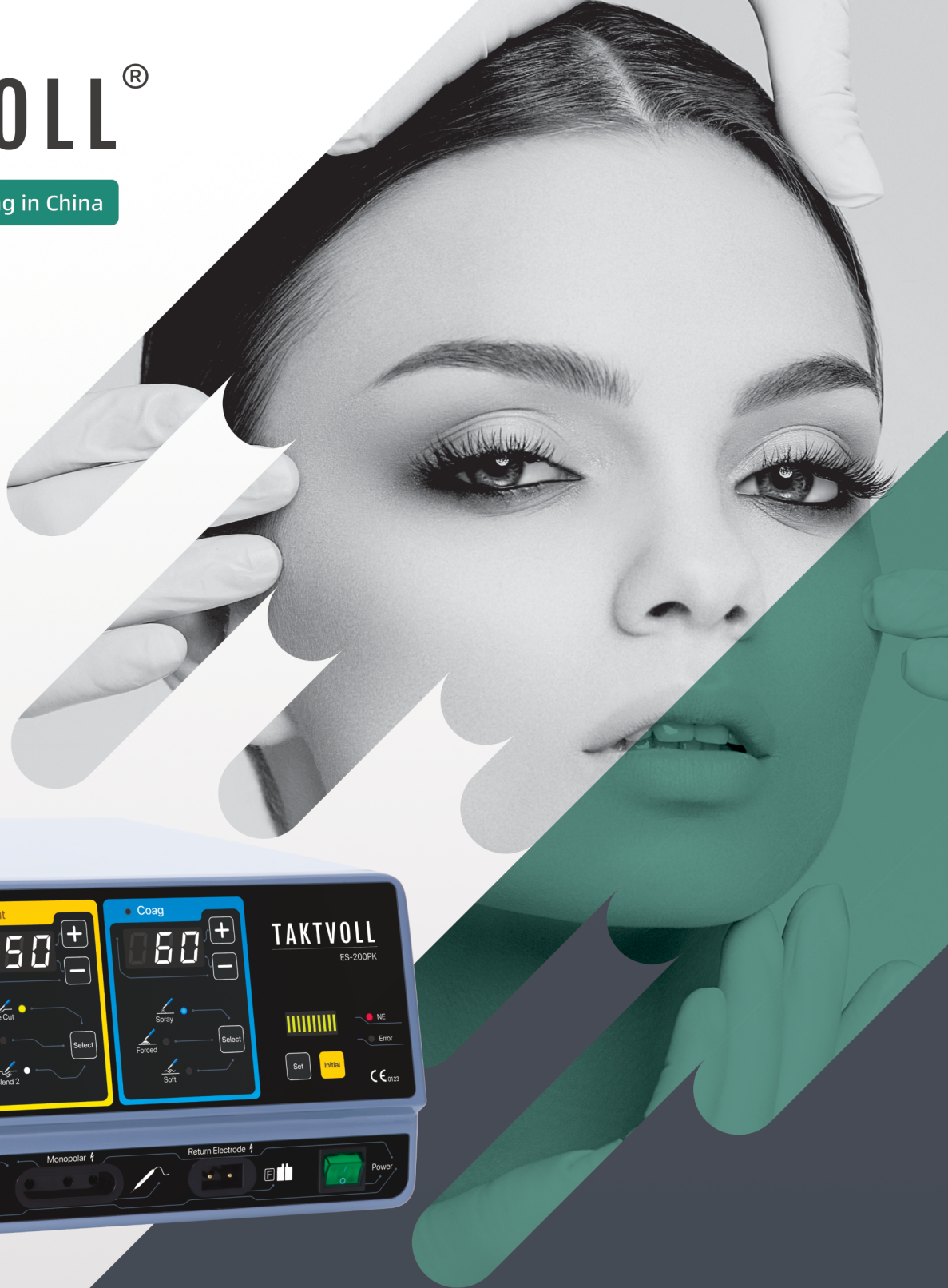

## Electrosurgical Generator Multifunctional Unit

ES-200PK

Taktvoll ES-200PK is a Multifunctional electrocautery generator with a wide range of application departments and a very high cost performance.

# Electrosurgical Generator

## ES-200PK Multifunctional Unit

The surgeon brings convenience and reduces surgical damage, and is particularly suitable for surgical applications such as general surgery, orthopedic surgery, gynecological surgery, ENT surgery, neurosurgery, skin plastic surgery, oral and maxillofacial surgery.



## ES-200PK Feature

### 3 Cut Mode

#### PURE

Provides a standard Cutting effect.

#### BLEND 1

Provides a combination of 70% Cutting and 30% Coagulation

#### BLEND 2

Provides a combination of 50% Cutting and 50% Coagulation

### 3 Coag Mode

#### SPRAY

Designed to provide hemostasis in a broader area with less depth of tissue damage when compared to desiccation.

#### FORCED

Fast and effective standard coagulation with moderate thermal penetration. Slight carbonization effects.

#### SOFT

Gentle coagulation with deep penetration, without carbonization, resulting in minimum adhesion of the electrode. Supported by the power control.

### 2 Bipolar Mode

#### Standard

Bipolar coagulation which gently coagulates tissue for minimal tissue damage reducing carbonization and sticking.

#### Precise

For microsurgical coagulation in the lower power range. Precise power settings and effects.

**TAKTVOLL®**

Beijing Taktvoll Technology Co., Ltd.

Address: Room 10a-302, No. 13, South Fourth Street, Jingsheng, Jinqiao Science And Technology Industry Base, Tongzhou District, Beijing, China.

Tel: +86(10) 56495579 / 56495580 Web: [www.taktvollmed.com](http://www.taktvollmed.com) E-mail: [Support@taktvoll.com](mailto:Support@taktvoll.com)