

TAKTVOLL®

Intelligent Manufacturing in China



PES-100PM


Electrosurgical System Portable LCD Touch Screen

Designed for precision, safety, and reliability, the PES-100PM Plus enables seamless monopolar and bipolar operation with dual footswitch interfaces for an efficient surgical workflow.



13485:2016

Portable 7-Inch LCD Touch Screen Electrosurgical System

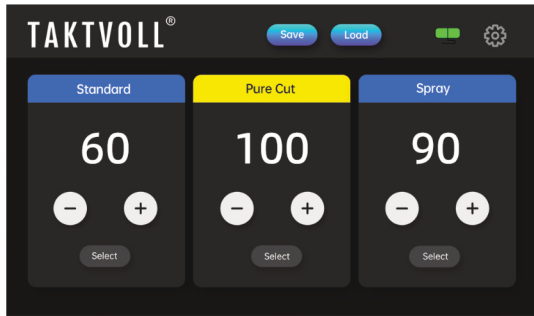

- Glass full version design, which is extremely tense and beautiful
- Color touch screen operation panel, flexible and easy to operate
- 7 Working modes
- Electrosurgical Pens and Foot Switch Control
- Recently Mode, Power and Other Parameters
- Quick Adjustment of Power
- Cut and Coag in an Intermittent Manner

WEIGHT AND DIMENSIONS	
Height	200mm
Width	210mm
Depth	175mm
Weight	2.0kg



PES-100PM Feature

3 Cut Mode		2 Coag Mode		2 Bipolar	
PURE	Provides a standard Cutting effect.	SPRAY	Designed to provide hemostasis in a broader area with less depth of tissue damage when compared to	Standard	Bipolar coagulation which gently coagulates tissue for minimal tissue damage reducing carbonization and sticking.
BLEND 1	Provides a combination of 70% Cutting and 30% Coagulation	FORCED	Fast and effective standard coagulation with moderate thermal penetration. Slight carbonization effects.	Precise	For microsurgical coagulation in the lower power range. Precise power settings and effects.
BLEND 2	Provides a combination of 50% Cutting and 50% Coagulation				



TAKTVOLL®

Beijing Taktvoll Technology Co., Ltd.

Address: Room 302, 3rd Floor, Building 10, No. 13, Jingsheng South 4th Street, Beijing Economic and Technological Development Zone, Beijing, China.

Tel: +86(10) 56495579 / 56495580 Web: www.taktvollmed.com E-mail: Support@taktvoll.com.